### Introduction

Human body generates very weak electricity that is not easily recognized. However it very sensitively responds to external electric stimulation. Electric treatment is to use this character to treat illnesses, pains and many other problems. Accordingly, if the low frequency waves are made to flow properly through the adjustment of strength and speed, they generate very rhythmical and smooth stimulating effects.

### Exercise Guidance

Within the first week you should give your muscles time to adjust to the new work-out. Use the unit only 2-3 a week for a maximum 10-15 minutes per muscle group. Do not increase the PULSE INTENSITY during the first week. After that you can use the SlenCare Butterfly for a maximum of 30 minutes per day. We recommended regular physical exercise in conjunction with treatment.

### How does the SlenCare massager work?

The Butterfly generates small, gentle pulses that delivered inside the pad placed on your skin. These safe pulses activates your underlying motor nerves, resulting in muscle contractions. This is the same technology as used in medicine, physiotherapy and sport. The massager will help to achieve total help, improve figure problems, especially flabby underarms, hips, sagging muscles to their firm. Regular exercise and healthy diet may improve you cardio-respiratory system. Naturally, results will vary from individual to individual, depending on a lot of factors, such as the individual's lifestyle, physiology and etc.

### Features

Relieve the pain in shoulders, back, waist and other areas of the body.  
Soothe anaesthetization and reduce fatigue.  
Help blood circulation.  
Beautify the skin and restore skin elasticity.  
Relieve muscle pain and increase fiber contraction.

### Functions

6 mode options from A to F with LED display.  
10 level strength.  
Automatically turn off after 12 min. operation.  
4 LED for indication

### Precautions

Attention! Read the manual thoroughly before use the product.  
The SlenCare Butterfly is not a medical device and recommended for home use. The massager extremely safe and easy to use, however it does not take any responsibility for any harm caused due to misuse of the product. This device should not be used in the following cases:  
tuberculosis,  
benign and malignant tumors,  
inflammations of the skin, phlebitis and thrombosis,  
open and fresh wounds,  
heart diseases  
broken skin or varicose veins.  
Pregnant women should avoid using the device.

*In case of any problems with health, please consult your doctor before using the massager.*


### Warnings

Please stop using it and consult your physician when unpleasant feeling occurs in the use.  
Please do not apply the device in your mouth, over the head, heart or on any damaged skin areas.  
Please shift all controls to the off positions before use.  
Please do not put the device on necklaces, rings or any other conductive metals.  
Please do not stick electrodes to each other.  
Please do not touch the massage pads with hands.  
Please do not use while driving.  
Please do not use while bathing or in a shower.  
Supervision is necessary when the device is used by children or disabled persons.

### Operating Instructions

#### Startup

1. Ensure the connection.
2. If the device is used for the first time, remove the battery case cover, insert the battery into the case with the "+" polarity facing up.
3. Clean the area where the massage pads will be attached.

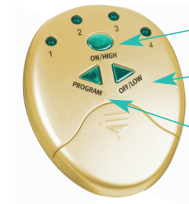- ON/HIGH** Button: Turns the unit on and increases the output intensity
- OFF/LOW** Button: Turns the unit off and decreases the output intensity. Each time this button pressed, the intensity decreases by one step down to level 0.
- PROGRAM** Button: This button is used to select one of the 6 pre-programmed exercise modes. The setting is indicated by color and location of the 4 LED lights on the unit.

4. Remove pad protector and keep it for later use.
5. Attach the device to the area to be massaged.
6. Press the ON/Hi button to turn on the power.

### Change the Modes

1. Press "Program" button once, the mode is switched from current one to the next mode sequentially. The setting is indicated by color and location of the 4 LED lights on the unit.
2. Mode change won't affect the stand-by status, working status and the strength as well.

### Select Mode

Mode A	LED 1 on	Tapping
Mode B	LED 2 on	Rhythmical Kneading
Mode C	LED 1 & 2 on	Deep Kneading
Mode D	LED 3 on	Pressing
Mode E	LED 4 on	Slow Stroke
Mode F	LED 3 & 4 on	Continuous Pressing

### Basic Operations

1. Press the "ON/HIGH" button again, the stand-by status is switched to working status.

### Change the intensity

1. If you want to increase the intensity level, use ON/HIGH button. Each press brings the strength to next higher level until level 10 is reached.
2. If you want to decrease the intensity level, use OFF/LOW button. Each press on "OFF/LOW" button brings to the next lower level until stand-by status is reached. One more press turns the device off.

### Switching-Off

1. To switch off the device, press ON/LOW until the device turns off.
2. When the device is turned off on a certain mode, the mode is saved for future use.

### Cleaning and Maintenance

Massager is designed for household use and is virtually maintenance free. No lubricants are needed. To clean use a damp cloth only. Do not use any strong chemical cleaners. Please put the massager back to cover after use. Store the massager in a safe place away from sharp objects. The device must be kept out of reach of children!

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SAVE THESE INSTRUCTIONS

For further information please contact:

SLENCARE LTD  
Suite 38-40  
London Fruit Exchange  
Brushfield Street  
London E1 6EU  
United Kingdom

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